OPTIKAM PRO 3LT, PRO 5LT, PRO 8LT

INSTRUCTION MANUAL


| Model |
|-----------|
| 4083.11LT |
| 4083.12LT |
| 4083.18LT |

Version: 1
Issued: 08, 11, 2016

M-005

Micrometer slide for software calibration, 1mm/10um, 10mm/100um.

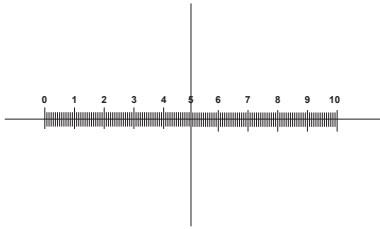1 DIV=0.1mm


For stereo microscopes calibration

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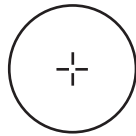1 DIV=0.01mm


For biological microscopes calibration

For biological microscopes calibration

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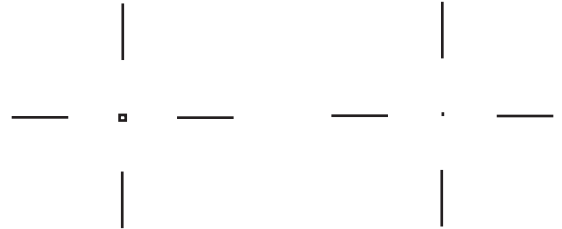Ø 0.15mm

Ø 0.07mm

For checking the calibration done by drawing a line or a circle

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Installation of camera drivers and software

Operating system requirements: Windows XP, Windows Vista, Windows 7, Windows 8 or Windows 10. If your PC is Windows XP : before the installation of your Camera please check whether it has Service Pack2 or not, should not please update your XP from the Microsoft web site. If your PC has got a USB less than 2.0 there may be some problems in capturing image from the camera due to the transmission not as fast enough.

1. please insert your CD and choose the folder \Optika Vision Pro\
please select the \camera drivers\ folder.
run vcredist_x32_bit.exe or vcredist_x64_bit.exe accordingly to your Windows PC characteristics.
run optikam_pro_plus_x32_bit.exe or optikam_pro_plus_x64_bit.exe accordingly to your Windows PC characteristics.
now please select the \ovp_4_5\ folder and run the setup.exe
2. The installer may show warnings about missing Windows certification, ignore such warnings and select "Continue". Note: do NOT cancel the installation.
3. Then Optika Vision Pro is installed.
4. At the end of the installation process restart the PC if requested.
You will find on your Desktop some icons which call from the just installed Optika software, at the same time all these icons will appear in the Start-programs-menu.
5. Plug-in the Camera in your PC, the camera will be identified and Windows will recognize the driver which has just been downloaded from the CD.
Note: it is not necessary to REinsert the CD again.
Then ignore any message about Windows compatibility and select "Continue" .
Note: do NOT cancel the installation.
6. At the end of the identification process, plug-off your Camera and try to plug-in it again, for checking if the Driver was installed well no Windows error messages should appear.
Should any error message appear on your PC, you will find a \Troubleshooting\ folder on your CD which shows how to recover the PC from the problems which Windows could give.

Some notes

You should know that your Windows PC needs to install a Driver before capturing images from the Camera, when the CD installs the Driver it won't be ready to be used at once, in fact Windows needs to recognize the camera before, and then it will finally put the Driver in the right \driver-folder\.

Note that in any PC there are couples of USB ports, should you connect only the Camera in a couple of port is better because you'll get a faster communication.

Often a PC has some ports on the front and some on the back, well, the ports on the back are usually faster.

If you connect the Camera into another USB port Windows needs to recognize its Driver again even if it has already done this operation for the first USB port.

Generally, every time you use a peripheral device and after connect your Camera then windows needs to recognize the Driver again.

It may happen sometimes that Windows doesn't recognize a Camera using one of your PC USB port and no problem with another one.

To prevent any installation problem we have put on the CD all the software in a specific \folder\ which has the same name.

How to test your Camera

For checking if your Camera works well is not necessary put it into the Microscope, it's enough to connect into the USB port, run Optika Vision Pro using the icon on your Desktop and check whether you can see a NOT black image.

What is your Driver

For Optikam Pro3 Camera the driver is:

Optika 319cu

For Optikam Pro5 Camera the driver is:

Optika 519cu

For Optikam Pro8 Camera the driver is:

Optika 891cu

Using Optika Vision Pro you don't have to tell the program how is the the Driver because it will be automatically recognized.

Getting started with Optika Vision Pro, a very complete measuring and documentation software.

Run Optika Vision Pro by selecting "Start", "Programs", "Optika Microscopes".

You can see the live image from the Camera, and you can now choose the right CCD exposure time, the gain, the colour gain, the white balance or the Gamma, in order to see the image in many different modes !


By selecting on the "camera" icon you can fix the current image in order to manipulate it by the software.

As you can see there are many measuring options on the panel: make a calibration, take segments, areas, perimeter measures, angles, ecc..

About how they work you can take a look at the manual on the Help menu or at the same into the CD.

Equipment disposal

Art.13 Dlsg 25 July 2005 N°151. "According to directives 2002/95/EC, 2002/96/EC and 2003/108/EC relating to the reduction in the use of hazardous substances in electrical and electronic equipment and waste disposal."

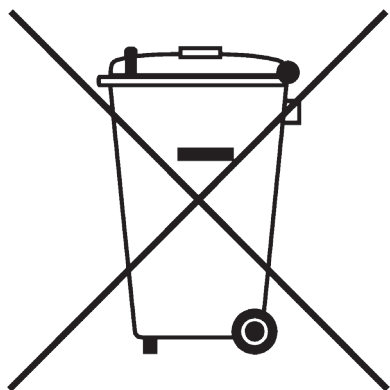The basket symbol on equipment or on its box indicates that the product at the end of its useful life should be collected separately from other waste.

The separate collection of this equipment at the end of its lifetime is organized and managed by the producer. The user will have to contact the manufacturer and follow the rules that he adopted for end-of-life equipment collection.

The collection of the equipment for recycling, treatment and environmentally compatible disposal, helps to prevent possible adverse effects on the environment and health and promotes reuse and/or recycling of materials of the equipment.

Improper disposal of the product involves the application of administrative penalties as provided by the laws in force.

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OPTIKA S.r.l.